



Conventional Laryngoscopes

No Compromise Professional Instruments



proactmedical.co.uk



“Supplying technology with no compromise
for over 20 years”

METAL^{MAX}TM

METAL^{MAX}TM
90

METAL^{MAX}TM
Combi

Metal Max™ 90 Disposable Laryngoscope Blades



- Conventional type fitting compliant to ISO 7376:2009
- The third generation of the original metal through the hook blade. The Metal Max™ 90 range was the world's first laryngoscope blade to feature an economical plastic block with metal through the critical hook on area. Where PROACT leads, others follow!
- Premium metal encased meLED™ light source offers a cool, super bright light and protects the LED from damage for gold standard dependability
- Correctly formed smooth atraumatic tip with no sharp edges
- Superb light output with no heat output from the meLED™
- Full length blade spatula and ET tube guide surface - no compromise
- Supplied clean with no sterile shelf life issue
- Two sprung ball bearings give a secure fit and option of re-use on the same patient plus allowing multiple testing for use on crash trolleys etc
- Supplied in new colour coded easy open sachets for faster access and ease of use
- Guaranteed to fit all compliant conventional handles
- All blade dimensions can be found at www.proactmedical.co.uk/bladedims

Part No.	Fitting	Size	Length (mm)
MMCV90MAC1	Conventional	MAC 1	75.0
MMCV90MAC2	Conventional	MAC 2	93.0
MMCV90MAC3	Conventional	MAC 3	110.0
MMCV90MAC4	Conventional	MAC 4	135.0
MMCV90MIL00	Conventional	MIL 00	43.5
MMCV90MIL0	Conventional	MIL 0	55.0
MMCV90MIL1	Conventional	MIL 1	80.5
MMCV90MIL2	Conventional	MIL 2	130.5



Metal Max™ 90 Disposable Lever Tip Laryngoscope Blades



- All the benefits of the Metal Max™ 90 blade as above - based on MAC blades
- Features a soft action Lever Tip function that assists with difficult intubation
- Provides sensitive controlled tip movement useful in a wide range of difficult intubation scenarios, such as:
 - Posterior tongue displacement
 - Prominent/unusual dental structures
 - Neck/spinal injuries
 - Recessive mandible
 - Forward displacement of the larynx

Part No.	Fitting	Size	Length (mm)
MMCV90LT2	Conventional	Lever Tip 2	93.0
MMCV90LT3	Conventional	Lever Tip 3	110.0
MMCV90LT4	Conventional	Lever Tip 4	135.0



Metal Max™ 100 Disposable / Speciality Laryngoscope Blades



- Conventional type fitting compliant to ISO 7376:2009
- Top of the range 100% stainless steel blade with same blade construction as a full reusable the original and best full metal disposable Laryngoscope Blade!
- Premium metal encased meLED™ light source offers a cool, super bright light and protects the meLED™ from damage for gold standard dependability
- Correctly formed smooth atraumatic tip with no sharp edges
- meLED light source gives out very little heat adding to patient comfort and safety
- Full length blade spatula and ET tube guide surface - no compromise
- Supplied clean with no sterile shelf life issue
- Triple sprung ball bearings give a secure fit
- Supplied in easy open sachets for faster access and ease of use
- Guaranteed to fit all compliant conventional handles
- All blade dimensions can be found at www.proactmedical.co.uk/bladedims



Wisconsin



Robertshaw



Cardiff Pro



Seward



Polio



Magill

Part No.	Fitting	Size	Length (mm)
HMM100MAC0	Conventional	MAC 0	55.0
HMM100MAC1	Conventional	MAC 1	75.0
HMM100MAC2	Conventional	MAC 2	96.0
HMM100MAC3	Conventional	MAC 3	113.0
HMM100MAC4	Conventional	MAC 4	137.0
HMM100MAC5	Conventional	MAC 5	153.0
HMM100MIL00	Conventional	MIL 00	47.5
HMM100MIL0	Conventional	MIL 0	58.5
HMM100MIL1	Conventional	MIL 1	82.5
HMM100MIL2	Conventional	MIL 2	134.5
HMM100MIL3	Conventional	MIL 3	175.0
HMM100MIL4	Conventional	MIL 4	185.0

Part No.	Fitting	Size	Length (mm)
HMM100WISC0	Conventional	Wisconsin 0	71.0
HMM100WISC1	Conventional	Wisconsin 1	90.0
HMM100WISC2	Conventional	Wisconsin 2	100.0

Part No.	Fitting	Size	Length (mm)
HMM100RS0	Conventional	Robertshaw 0	70.0
HMM100RS1	Conventional	Robertshaw 1	90.0

Part No.	Fitting	Size	Length (mm)
HMMLT100MAC3	Conventional	Lever Tip 3	108.0
HMMLT100MAC4	Conventional	Lever Tip 4	133.0

Part No.	Fitting	Size	Length (mm)
HMM100CP00	Conventional	Cardiff Pro 00	57.0
HMM100CP0	Conventional	Cardiff Pro 0	73.0
HMM100CP1	Conventional	Cardiff Pro 1	79.0

Part No.	Fitting	Size	Length (mm)
HMM100MAGILL0	Conventional	Magill 0	74.0
HMM100MAGILL1	Conventional	Magill 1	83.0

Part No.	Fitting	Size	Length (mm)
HMM100SEW1	Conventional	Seward 1	83.0
HMM100SEW2	Conventional	Seward 2	108.0

Part No.	Fitting	Size	Length (mm)
HMM100POL3	Conventional	Polio 3	118.0
HMM100POL4	Conventional	Polio 4	138.0

Metal Max™ Disposable Laryngoscope Handles



- Professional strength full 100% metal handle with grooves for good grip
- Conventional type fitting compliant to ISO 7376:2009
- Requires standard batteries to be fitted i.e x2 C size Adult or x2 AA Paediatric
- Available in a variety of sizes from adult, paediatric to stubby
- Supplied clean with no sterile shelf life issue
- Supplied in new colour coded easy open sachets for faster access and ease of use
- Eliminates the capital, sterilising and running costs of reusable handles
- Guaranteed to fit all compliant conventional blades

Part No.	Fitting	Size
MMCVH001	Conventional	Adult
MMCVH002	Conventional	Paediatric
See below	Conventional	Stubby





Metal Max™+ Disposable Laryngoscope Handles with Preloaded Batteries



- Professional strength full 100% metal handle with grooves for good grip
- As above, but with fitted 6V batteries (4x 1.5V in clear shrink wrap to make one cell)
- Fast and easy battery removal at end of life plus Eco friendly battery technology ensures ease of disposal and environmental efficiency
- Designed to be used only with Metal Max™ blades the 6V Metal Max+ handles offer an approx 33% increased light output ensuring optimum illumination
- Long battery life - over 5 hours continuous illumination. 3 years shelf life
- 31% lighter than the equivalent handle used with C size batteries

WARNING Only for use with Metal Max™ Blades

Part No.	Fitting	Size
MMCVHB01	Conventional w/batteries	Adult
MMCVHB02	Conventional w/batteries	Paediatric
MMCVHB03	Conventional w/batteries	Stubby

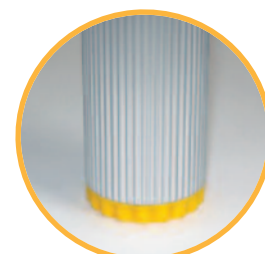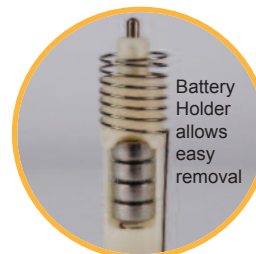


4x 1.5V in clear shrink wrap to make one easy to remove cell

Hydra Autoclavable Waterproof Reusable Laryngoscope Handles

- Hydra handles are waterproof and fully autoclavable (remove batteries first) guaranteed for up to 3000 cycles at up to 138°C
- Washer / disinfectant compatible (see user instructions at www.proactmedical.co.uk/support)
- Textured finish for a secure grip
- Lifetime Product Warranty

Part No.	Fitting	Size
HYCVH001	Conventional Waterproof	Adult
HYCVH002	Conventional Waterproof	Paediatric
HYCVH003	Conventional Waterproof	Stubby

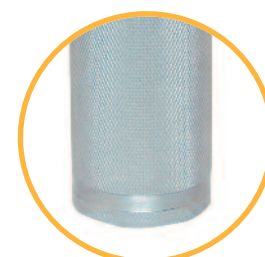

Metal Max™+ Disposable Combi™ Laryngoscope Sets - Prefitted 6V Batteries



Where PROACT leads, others now follow! Invented by PROACT, the Combi™ Laryngoscope concept provides an all in one Laryngoscope set with all the same features as the Metal Max 90 meLED Laryngoscope but with a pre-fitted PERMANENTLY attached blade. This critically, allows a Combi Set to be tested within its bag and eliminates any possibility of blades not fitting the handle or, misplacement of blade or handle. PROACT Combi™ Sets offer unparalleled reliability and performance and are ready to use in seconds.

Now also available in Lever Tip configuration, the Combi™ Laryngoscope range provides an all-in-one solution, including handle, 6V preloaded batteries and complete with a colour coded test-in-bag protective sheath. The simplest, fastest and most dependable solution available - try for yourself.

- **Ready for use in less than 5 seconds!**
- **Essential test-in-bag capability for repeated testing and real peace of mind**
- Non-removable Metal Max™ 90 blade permanently attached to a lightweight super strong aluminium full metal handle removes hassle and compatibility issues
- Super strong Metal Max™ blade with metal through the hook for maximum strength
- Full length blade spatula and ET tube guide surface - no compromise
- 31% lighter than the equivalent set when fitted with normal batteries
- Long battery life - over 5 hours continuous illumination. 3 years battery shelf life from the date of manufacture
- Fast and easy battery removal at end of life
- Eco friendly alkaline battery technology ensures ease of disposal and environmental efficiency
- Supplied in easy open colour coded packaging with no sterile shelf life issue
- All blade dimensions can be found at www.proactmedical.co.uk/bladedims



Part No.	Fitting	Size	Length (mm)
MMCB90MAC2	Combi™	MAC 2	93.0
MMCB90MAC3	Combi™	MAC 3	110.0
MMCB90MAC4	Combi™	MAC 4	135.0
MMCB90MIL00	Combi™	MIL 00	43.5
MMCB90MIL0	Combi™	MIL 0	55.0
MMCB90MIL1	Combi™	MIL 1	80.5

Part No.	Fitting	Size	Length (mm)
MMCB90LT2	Combi™	Lever Tip 2	93.0
MMCB90LT3	Combi™	Lever Tip 3	110.0
MMCB90LT4	Combi™	Lever Tip 4	135.0



ESSENTIAL TEST IN BAG PACKAGING

