



Gaumard[®]
Simulators for Health Care Education



HAL[®] S315.600

CPR + D Skills Trainer with OMNI[®] 2

- Wireless control with touchscreen OMNI[®] 2 from up to 30 feet
- Programmable spontaneous breathing and pulse
- Illustrates 20+ EKG rhythms on real monitors
- Deliver up to 360 joules of energy with real equipment
- eCPR[™] - CPR effectiveness monitor and trainer


HAL® S315.600 | CPR+D Skills Trainer with OMNI® 2

EASY-TO-USE, DURABLE, & PORTABLE HANDS-ON BASIC LIFE SUPPORT TRAINING

Whether for outreach or professional development, the new HAL S315.600 is a great addition to every basic life support algorithm training curriculum. Train learners on assessing breathing, performing quality CPR, delivering a shock, and recognizing return of spontaneous circulation. HAL offers a portable, all-in-one resuscitation training solution for every training need.



Give real-time feedback on CPR quality, or use coaching mode to learn

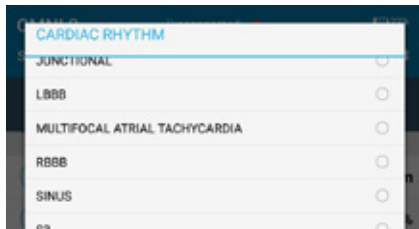

Illustrate shockable and non-shockable EKG with 20+ rhythms



Perform nasal and oral intubation with standard, real adjuncts




Visualize realistic chest rise with BVM ventilation




Identify proper chest rise and pulses without add-on items



Monitor and cardiovert rhythms using a real AED/defibrillator

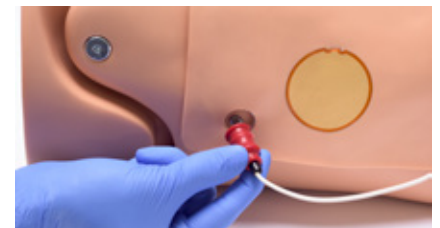


Defibrillate directly on patches or pad-free posts numerous times

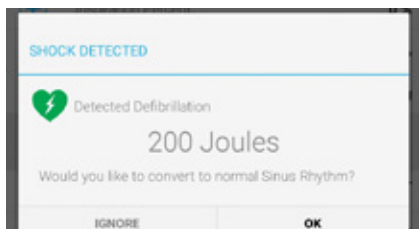

See immediate feedback on actions and select the responses



Check for return of spontaneous circulation

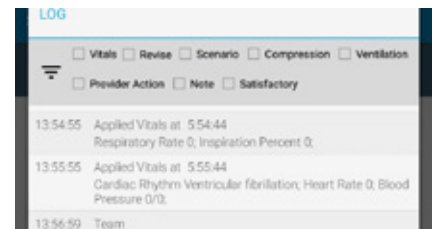

Automatically detect and record events for easy debriefing



DEBRIEF, LEARN, & IMPROVE

Save, print, and debrief performance reports to remediate, then repeat complex scenarios to prepare for real world emergencies.

RESULTS		
Total Session Time	00:07:04	
C/V Cycles Performed	29	
COMPRESSION RATE	COMPRESSION DEPTH	AVERAGE RECOIL
105	5.1	0.2
Target Compression Rate 100 CPM	Target Compression Depth 5.0 cm	Target recoil 0 in.
AVERAGE VENTILATION RATE	AVERAGE VENTILATION PIP	AVERAGE VENTILATION DURATION
10	25	1.3
Target Ventilation Rate 12 VPM	Target Ventilation PIP 17 cmH ₂ O	Target Ventilation Duration 1 sec
CLOSE		SAVE REPORT

AIRWAY AND BREATHING

- Practice nasal and oral intubation (ETT, Supraglottic Airway, King LT®) with realistic tongue, vocal cords, and trachea
- Illustrate head tilt/chin lift and jaw thrust
- Gastric distension with esophageal intubation and/or excessive ventilation
- Right mainstem intubation presents unilateral chest rise
- Visible chest rise with bag valve mask ventilation
- Spontaneous breathing (Automatic)
- Ventilation detection sensors

CARDIAC & CIRCULATION

- Anatomically correct landmarks for proper CPR hand placement
- Realistic chest cavity resistance and recoil
- Chest compression depth sensor
- Real-time CPR feedback
- Effective chest compressions generate palpable carotid pulses
- Monitor and cardiovert rhythms using a real AED/defibrillator
- Illustrate 20+ shockable and non-shockable EKG rhythms
- Deliver up to 360 joules of real energy to skin patches or snap connectors
- Palpable carotid pulses with variable rate synchronized with EKG

PERFORMANCE MONITORING

- Real-time CPR quality feedback
 - Compression depth and rate
 - Ventilation rate
 - Excessive ventilation
 - No-flow time
 - CPR cycles
- CPR metronome: audible tones help guide correct compression and ventilation rate and ratio
- Compliant with 2015 adult resuscitation guidelines
- Built-in resuscitation algorithm checklists for tracking individual and team actions

OMNI 2 FEATURES

- Built-in wireless connectivity with range up to 30 ft.
- Wireless interface with 8" touchscreen display
- Optimized vital sign controls for on-the-fly operation
- Built-in library with 20+ cardiac rhythms and options
- Compatibility with optional virtual vital sign to display 10+ parameters for delivering post-cardiac arrest care
- Alternate virtual shock function capability
- Session log records provider actions, vital signs, CPR metrics, and comments

DEBRIEFING TOOLS

- eCPR™ performance report provides averages for each CPR metric and cycle
- Save, email, and print eCPR performance reports for debriefing and archiving

CPR+D SKILLS TRAINER WITH OMNI® 2

S315.600

Available skin tones



PACKAGE CONTENTS

- HAL S315.600 Skills Trainer
- OMNI 2 Wireless Control Interface
- Carrying Case
- Operating Guide
- A/C Adapter

MODIFIED PHILIPS® DEFIBRILLATOR CABLE

S315.600.126

MODIFIED PHYSIO-CONTROL® DEFIBRILLATOR CABLE

S315.600.127

MODIFIED ZOLL® DEFIBRILLATOR CABLE

S315.600.128

REQUEST A QUOTE

14700 SW 136th Street
Miami, Florida 33196-5691

P.O. Box 140098
Coral Gables, Florida 33114-0098

Toll Free USA
Call 8:00 a.m. - 7:30 p.m. ET
Monday - Friday
800.882.6655

Worldwide
305.971.3790

Fax
305.252.0755
24 hours

Web Site
www.gaumard.com

Sales / Customer Service
sales@gaumard.com

ABOUT GAUMARD®

We are your direct source for health care education needs. Gaumard® is family owned and operated.

PLACE YOUR ORDER

By telephone, fax, or on-line
at www.gaumard.com

WE ACCEPT

VISA, MasterCard, American Express

WARRANTY

Gaumard products are covered by a one-year limited warranty. Terms and conditions apply. Please visit www.gaumard.com for details.

EXTENDED WARRANTY & TRAINING

Extended warranty plans and training are available.

TECHNICAL SUPPORT

Available 8:00 a.m. to 7:30 p.m. ET weekdays, or email us at support@gaumard.com
Phone: 786-478-3838

ON-LINE SUPPORT 24/7

Download articles, instructions, teaching tips. View instructional videos on demand.

REPAIR

Gaumard® offers repair services and parts. For more information visit www.gaumard.com

GAUMARD® LISTENS

Customers are a terrific source for ideas! If you need something special, give us a call or email us at sales@gaumard.com



Gaumard®
Simulators for Health Care Education