

Delivery bed

# Nascentia

The delivery bed that is comfortable,  
easy to use and ergonomic

Comfortable and ergonomic for the mother, easy to use by midwives,  
doctors and nursing staff. That characterizes the Nascentia delivery bed.



Overhead trapeze

Leg rests

Pelvic section (0-15 degrees)

Removable headboard

Electronic service panel

Firm grip

Remote control

Can be raised and lowered (60-100 cm)

Easy 360-degree turning radius

\*Electric drive system

\* Electric drive system is available as an accessory to facilitate transport.

# Delivery bed Nascentia

The Nascentia delivery bed is a comfortable, fully-equipped maternity and delivery bed, which makes the bed the perfect choice for the delivery room. The expecting, birthing and recovering woman can all make use of a single bed. The Nascentia is user-friendly, equipped with every essential feature and designed with all labor positions in mind.

By removing the leg section, the partus stool, the bed can be transformed into a delivery chair in one simple movement. This section can be used for various labor positions or as a chair for the doctor or obstetrician.



## Ergonomics

The Nascentia delivery bed is an ergonomic bed, with a focus on comfort and the prevention of back injuries, both for the mother, as well as for the nursing staff and doctors. The delivery bed can easily be handled by one person, even when it is occupied by a patient, due to the 360-degree turning radius.

# Features

The Nascentia delivery bed is characterized by its usability. Regulating the height, adjusting the footrests: everything is easy to operate. The Nascentia delivery bed has the following features:



## Naturally comfortable

- The foam mattress is specially designed for the delivery bed and offers the mother during and after the labor the desired comfort. The mattress has a density of up to 50 kg/m<sup>3</sup>.

## Real user-friendliness

- Accessories and features are easy to use and operate.
- The Nascentia has an integrated battery and can be operated without connection to the electrical grid.
- Due to clever wheel positioning, the bed has a 360-degree turning radius, which allows for easy control by one nurse.

## **Clean environment**

- No separate elements or options need to be installed. Every feature and accessor has been integrated into the bed's design.
- Leg supports are integrated under the bed.

## **Adjustable**

- Both the bed as a whole, as well as the leg section separately can be adjusted in height.

## **CPR-function**

- In case of resuscitation the CPR handle allows for a quick but gas spring controlled smooth flattening of the backrest.

## **Multiple sections**

- The various sections (the headrest, pelvic section and partus stool) can all be adjusted individually.
- The Nascentia delivery bed thus guarantees optimal positioning in every phase of the delivery.

## **Removable Partus stool**

- The Partus stool can be uncoupled from the bed and used as a chair for the doctor or nurse. It can be controlled manually when used separately, or automatically when connected to rest of the bed.

## **Safety**

- Due to the carefully designed columns, the bed offers high stability and support for the patient and the nursing staff.

## Labor positions



1. Throne position



2. Side Lying



3. Lithotomy



4. Kneeling birth position



5. Trendelenburg position



6. The sitting position





## Hygiene

The prevention of infections is a major health care issue. During the design process TechniCare naturally made sure to create a design that is in accordance with the regulations for infection prevention.

One of the most prominent features of the delivery bed is its cleanability.

The foam mattress has a washable polyester cover with PU coating for aeration, but is impermeable for any fluids. It is in compliance with the European fire standards EN1021-1/2 and EN 597-1/2.

## Research and Development


The delivery bed is developed in cooperation with multiple hospitals and distinguished suppliers. **“The perfect bed did not exist yet. So we decided to let it be designed for us”** said Tiba Spaapen, nursing unit manager at the OLVG Hospital in Amsterdam. **“We needed a single bed that improved the conditions both for the mother and the nursing staff.”**

This became a central idea during the design process.

## Healing environment

The delivery bed is designed according to the Healing Environment Concept, which means creating a customized environment that improves the mother’s well-being. A calming environment has a stress relieving effect on both the patient and the family.



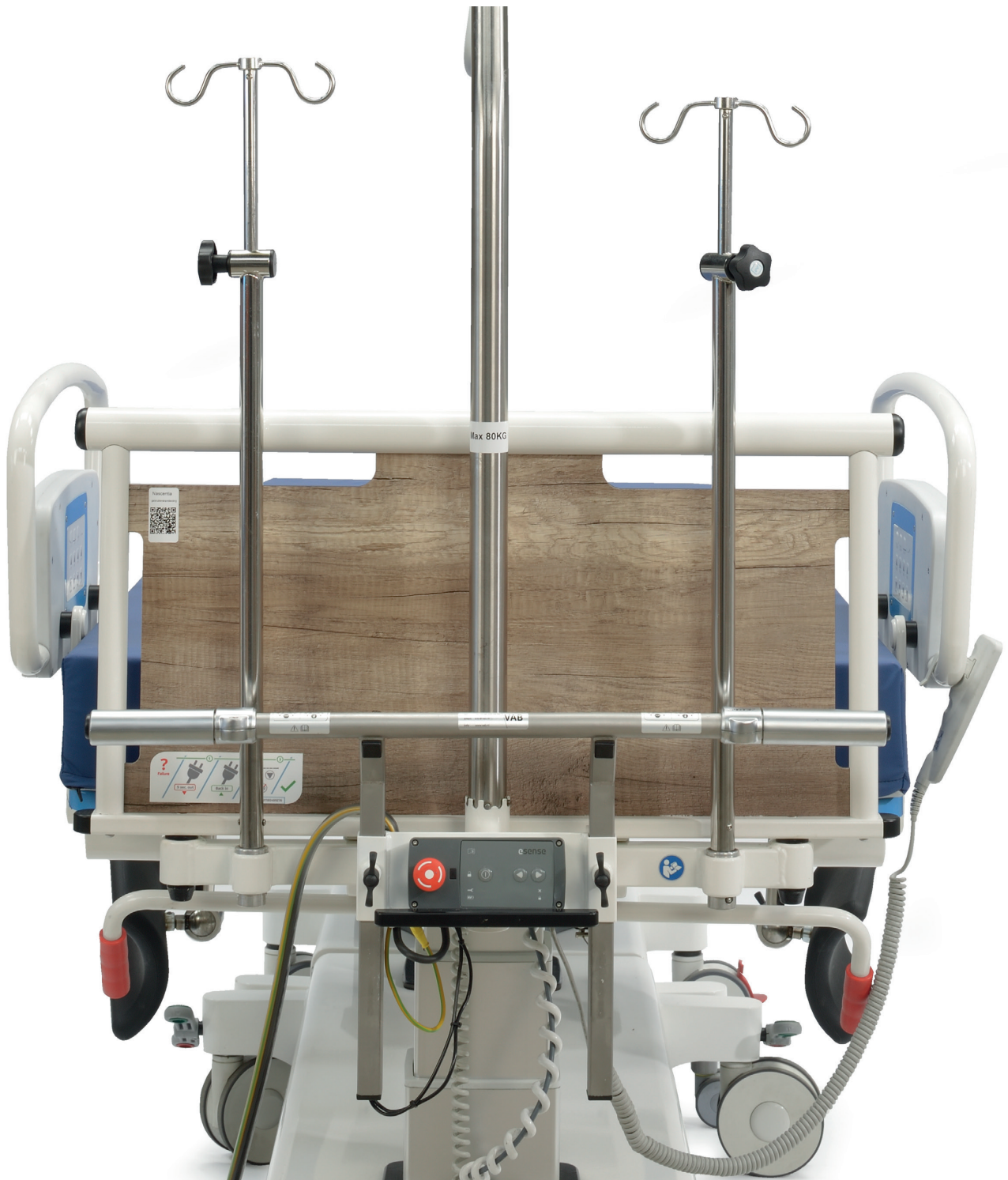


TechniCare delivers custom made accessories.  
Please ask your contact person at GM Medical  
for your desired accessoires.

## Colors



1.



2.



3.



4.



5.



6.



7.



8.



9.



10.



11.



12.



13.



14.



15.



16.



17.

1. Grey White

7. Grey Aluminium

13. Sapphire Blue

2. Pure White

8. Signal Yellow

14. Night Blue

3. Light Grey

9. Light Ivory

15. Turquoise (standard)

4. Antraciet Grey

10. Antique Pink

16. Pure Orange

5. Jet Black

11. Ruby Red

17. Turquoise Blue

6. White Aluminium

12. Yellow Green

## Accessories:

- Unique: Overhead adjustable Trapeze
- Under bed light
- Hand grip
- One extra IV pole
- Remote control
- V-cut mattress



## Technical specifications

Overall length:	235 cm
Overall width:	103 cm
Mattress length:	210 cm
Mattress width:	90 cm
Mattress thickness:	14 cm
Height adjustment:	60 - 100 cm
Adjustment backrest:	max 80 degrees
Pelvic section:	0-15 degrees
Lifting capacity - bed:	250 kg
Lifting capacity - seat:	200 kg (when detached)
Trendelenburg:	up to 15 degrees
Castors:	15 cm



**Maria Wallén**  
 Product specialist  
 Mobile +4673 588 63 01  
 mw@gm-medical.com